



PDXZ1(/M) Test Report

Date	2/20/2024	Production Order Number	W437473
Serial Number	43747301	Assembled By	LQ_Jiang
Range Tested	4.5 mm	Techician	LQ_Jiang


Forward Pitch/urad Vs. Position/mm

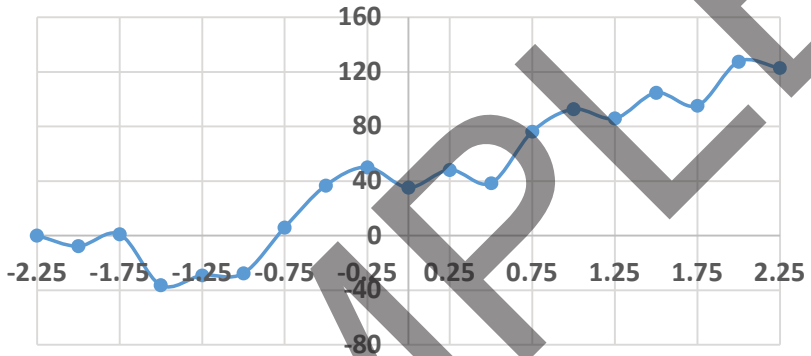


Backward Pitch/urad Vs. Position/mm

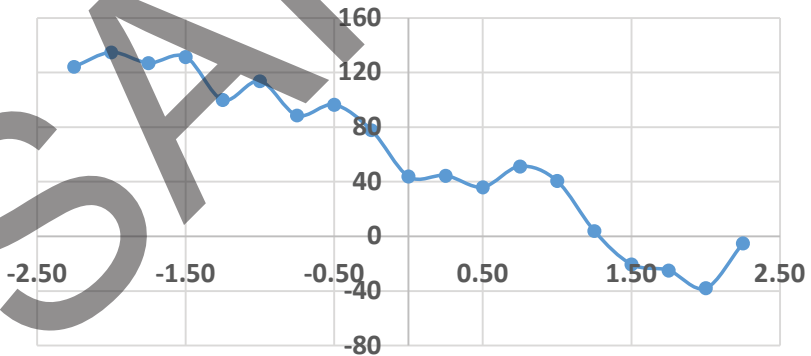

Forward Roll/ μ rad Vs. Position/mm

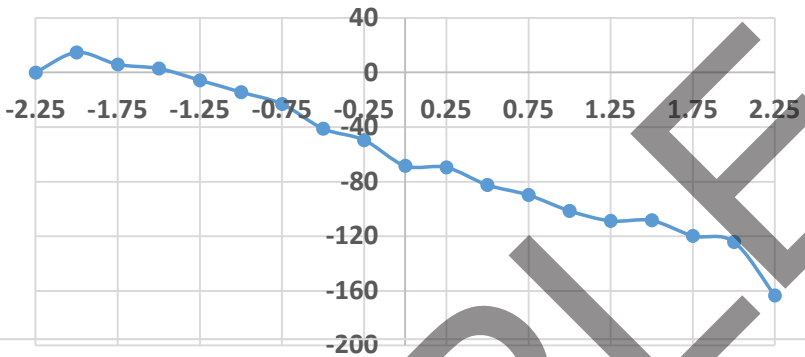

Backward Roll/ μ rad Vs. Position/mm

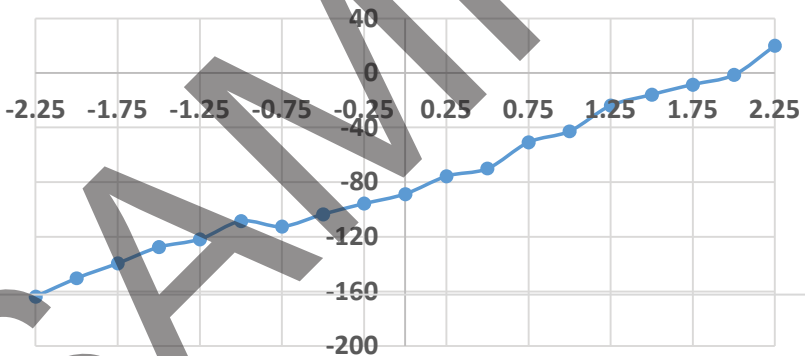

Items	Pitch/ μ Rad	Roll/ μ Rad
Forward	164	178
Backward	173	184
Specification	200	200